

说明书

材料：105哑粉

工艺：啤,风情7折

尺寸108x96mm

96

108

How to connect the App

Note: Do not use the device unless you have successfully connected the device to your mobile phone through the Android/iOS App and are streaming live video footage.

1

Download the APP

iOS: Search for "bebird" in App Store
Android: Search for "bebird" in Google Play




OR scan the QR code to download the app



2

Open the cover
Power on automatically

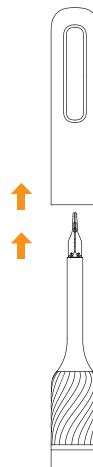


3

Open the APP
Connect device WiFi




4

Connect product WiFi
"bebird-Note5-xxxxxx"

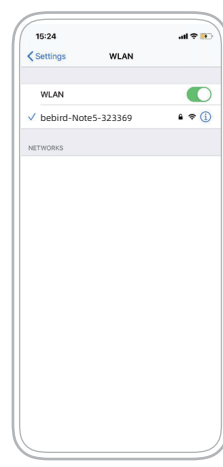


5

Back to bebird APP
Click on "GO" to use it

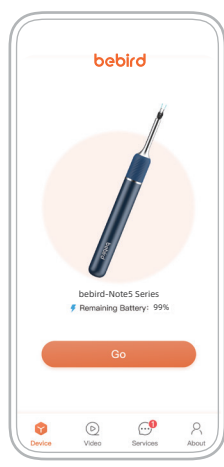

6

The Visual Screen

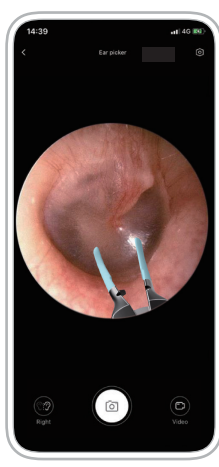


7

Press the button
to control tweezers





Product Overview

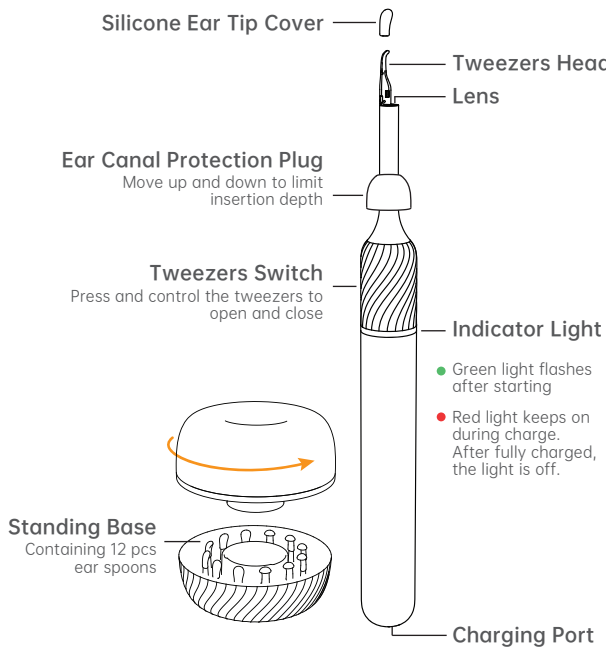

How to use the Tweezers


Hold the tweezer switch still with 2 fingers
and remove the ear spoon with the other hand.

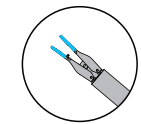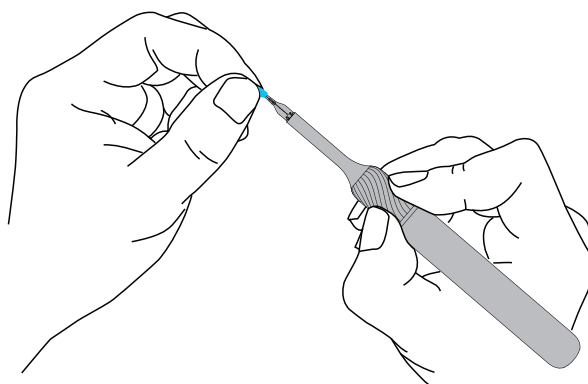

The tweezers open automatically
when you release your hand.

How to Replace Ear Spoon

Press the tweezer switch first so that the tweezers
are closed. Just put the ear spoon set on.

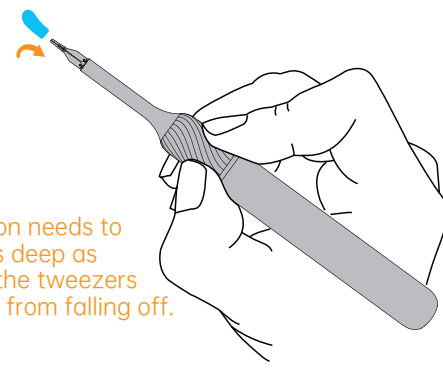

Note:
The ear spoon needs to
be placed as deep as
possible on the tweezers
to prevent it from falling off.

! If the ear spoon is damaged, please stop using
it immediately and contact us for replacement.

Accessories

1



Spatula ear scoop x 4
Multiple ear scoop covers,
used separately

2



Dome ear scoop x 4
Easily remove earwax in 360 degree

3




Sponge ear
spoon cover X4
Clean and comfortable,
remove earwax and clean oil ear

4




Ear Canal Protection
Plug x 1
A touch reminder to avoid
reaching too deep in the ear canal.

5

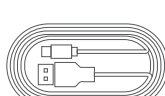

Charging Cable x 1

6

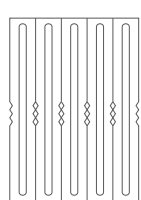

Alcohol swab x5
Cleaning the lens or spoon

7



Double-ended
cleaning brush x 1
Suede head to wipe the lens
Brush head to clean the ear tip

Troubleshooting

Problem	Solution
No response after startup	Please use the incidental charging cable to charge
The LED light of the camera is not on	Interrupt the connection, restart the product to connect again
The app exits quickly	Delete the existing app and scan QR to reinstall it
The indicator light blinks quickly	Please use the incidental charging cable to charge
Abnormal display in app	Interrupt the connection, restart the product to connect again
The Wifi connection was unstable and disconnected in a few seconds	Interrupt the connection, restart the product to connect again
The visual interface blur	Please wipe it with a double-ended cleaning brush or alcohol swab
Disorientation when in use	Please switch the left / right ear mode under the APP interface
If you still have problems with the above procedures, please contact our Customer Service Representatives directly for further assistance.	

Attention

- Please use the equipped alcohol swab to wipe the lens carefully when cleaning.
- The ear spoon cover must be installed properly, not loose before use to prevent it from falling off.
- Please observe the surrounding environment before use. Do not use it in an environment that is susceptible to impact.
- Prohibit children from using this product alone to avoid inadvertent injury.
- Please do not put the device inside any form of liquid.
- The battery must be charged once a month for effective use.
- The temperature will rise slightly when using.
- The product temperature will rise slightly when using (from bebird laboratory data, the maximum temperature of the whole machine does not exceed about 45 °C). Please feel free to use.

Warranty card

User information	
User name:	Telephone number:
Seller:	Shipping address:

Product parameter

Brand	bebird
Product Name	3 in 1 Smart Visual Ear Cleaner
Model	Note5 Pro
Pixel	10 Megapixels
Lens diameter	3.2 mm
Optimum focal length	13.5 mm
Depth of field	10 ~ 50 mm
Gyroscope	Gimballed Intelligent Gyroscope
Battery	230 mAh
Working time	About 90 minutes
Charging time	About 60 minutes
Image transmission rate	20 fps
Working temperature	-10 ~ 50°C
Input parameters	DC 5V ≈ 0.5A
Rod weight	92 g
Working frequency	2.4 Ghz
Network standard	IEEE 802.11b/g
Image sensor	CMOS



Hazardous substance

Part name	Hazardous material					
	(Pb)	(Hg)	(Cd)	(Cr(VI))	(PBB)	(PBDE)
Printed Circuit Components	×	○	○	○	○	○
Standing Base	○	○	○	○	○	○
Cover	○	○	○	○	○	○
Plastic rubber accessory	○	○	○	○	○	○

This form is formulated in accordance with SJ/T 11364
○ : represent the hazardous material ratio is below GB/T26572-2011
× : represent the hazardous material ratio is over GB/T26572-2011

bebird

3 in 1
Smart Visual
Ear Cleaner

Visible
Pickable
Clipable

⚠ Note: Do not use the charging equipment
with output voltage exceeding 5V